

Consort

MANUAL

HANDLEIDING

MODE D'EMPLOI

ANLEITUNG



Distributor:
Topac Inc
231 CJC Highway
Cohasset
MA 02025
Tel 781 740 8778
Fax 781 740 8779
www.topac.com
[email sales@topac.com](mailto:sales@topac.com)

C5010
C5020

Table of contents

| | |
|--|---|
| Introduction | 1 |
| Manufacturer | 1 |
| Warranty | 1 |
| Servicing..... | 1 |
| Batteries | 1 |
| AC adaptor | 1 |
| Keyboard | 1 |
| Set-up..... | 1 |
| GLP..... | 2 |
| Logging data | 2 |
| Recalling data | 2 |
| pH measurement..... | 2 |
| mV measurement | 2 |
| Inputs | 2 |
| Important | 2 |
| Conductivity measurement..... | 3 |
| TDS/Salinity measurement | 3 |
| Oxygen measurement | 3 |
| %O ₂ measurement | 3 |
| Temperature measurement | 3 |
| Maintenance of pH electrodes..... | 4 |
| Maintenance of metal electrodes | 4 |
| Maintenance of conductivity electrodes | 4 |
| Maintenance of oxygen electrodes..... | 4 |

Table de matières


| | |
|--|----|
| Préface | 9 |
| Fabricant | 9 |
| Garantie | 9 |
| Service..... | 9 |
| Batteries | 9 |
| Adaptateur secteur..... | 9 |
| Le clavier | 9 |
| Configurer..... | 9 |
| Branchement | 10 |
| Important | 10 |
| Mesure de pH..... | 10 |
| Mesure de mV | 10 |
| BPL | 10 |
| Mémoriser des données | 10 |
| Afficher des données..... | 10 |
| Mesure de conductivité | 11 |
| Mesure de TDS/salinité | 11 |
| Mesure d'oxygène | 11 |
| Mesure de %O ₂ | 11 |
| Mesure de température..... | 11 |
| Entretien des électrodes de pH | 12 |
| Entretien des électrodes métalliques | 12 |
| Entretien des électrodes de conductivité ... | 12 |
| Entretien des électrodes d'oxygène..... | 12 |

Inhoudstafel

| | |
|---|---|
| Inleiding | 5 |
| Fabrikant | 5 |
| Waarborg | 5 |
| Onderhoud | 5 |
| Batterijen | 5 |
| Netadapter | 5 |
| Instellen | 5 |
| Toetsenbord..... | 5 |
| Belangrijk | 6 |
| Gegevens oproepen..... | 6 |
| GLP | 6 |
| mV-meting | 6 |
| pH-meting..... | 6 |
| Gegevens opslaan..... | 6 |
| Ingangen..... | 6 |
| Geleidbaarheidsmeting | 7 |
| TDS-/Saliniteits-meting | 7 |
| Zuurstofmeting | 7 |
| %O ₂ -meting | 7 |
| Temperatuurmeting | 7 |
| Onderhoud van pH elektroden..... | 8 |
| Onderhoud van metaalelektroden | 8 |
| Onderhoud van geleidbaarheidscellen | 8 |
| Onderhoud van zuurstofelektroden..... | 8 |

Inhaltsverzeichnis

| | |
|---|----|
| Hersteller | 13 |
| Garantie | 13 |
| Service | 13 |
| Einführung | 13 |
| Batterien | 13 |
| Netzadapter | 13 |
| Bedienungsfeld | 13 |
| Einstellen..... | 13 |
| Eingänge..... | 14 |
| Wichtig | 14 |
| pH-Messung | 14 |
| mV-Messung | 14 |
| GLP | 14 |
| Daten speichern..... | 14 |
| Daten auflisten | 14 |
| LF-Messung | 15 |
| TDS-/Salinität-Messung | 15 |
| Sauerstoffmessung | 15 |
| %O ₂ -Messung..... | 15 |
| Temperaturmessung | 15 |
| Wartung der Metallelektroden | 16 |
| Wartung der pH-Elektroden | 16 |
| Wartung der Leitfähigkeitsmeßzellen | 16 |
| Wartung der Sauerstoffelektroden..... | 16 |

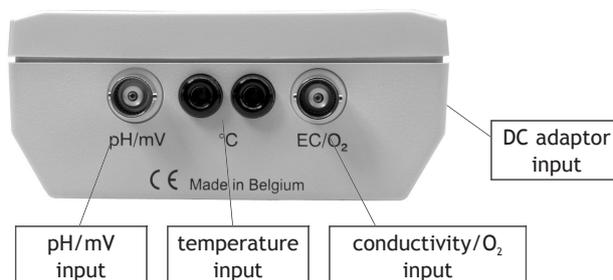

Introduction

This instrument is manufactured with the latest technology and needs no particular maintenance. **Consort** certifies that this instrument was thoroughly inspected and tested at the factory prior to shipment and found to meet all requirements defined by contract under which it is furnished. However, dimensions and other physical characteristics may differ.

The normal operating temperature should be between 0° and 40°C. Never use the instrument in a room with high humidity (>95 %) or at very low temperatures (condensation water!).

Manufacturer

Consort bvba

Parklaan 36

B2300 Turnhout

Belgium

Tel (+32)(14)41 12 79

Fax (+32)(14)42 91 79

E-mail: info@Consort.be

Warranty

This instrument (excluding all accessories) is warranted against defective material and workmanship for a period of thirty-six (36) months from the date of shipment ex factory. **Consort** will repair all defective equipment returned to it during the warranty period without charge, provided the equipment has been used under normal laboratory conditions and in accordance with the operating limitations and maintenance procedures in this instruction manual and when not having been subject to accident, alteration, misuse or abuse. A return authorisation must be obtained from **Consort** before returning any product for warranty repair on a freight prepaid basis!

Consort is not liable for consequential damages arising out of the use or handling of its products.

Servicing

In the event of this instrument being returned for servicing, the owner is requested to remove the power supply lead and **NOT** to send the following items unless they are suspect:

Manual

Cables

Accessories

If serious malfunctioning occurs, stop using the unit immediately and consult your local **Consort** dealer.

Batteries

The batteries can be replaced by opening the bottom compartment of the cabinet. Only use NiMH batteries!

AC adaptor

Connect the AC adaptor (9 V, 300 mA for 230 V~ or 120 V~) to the DC socket. **Avoid to use an AC adaptor without NiMH batteries being inserted in the instrument!** Do not hold the adaptor by wet hand!

Keyboard

MODE = Selects all modes or escapes from error traps, calibration procedures, etc. by returning to the original mode.

CAL = Starts or proceeds a calibration or a function.

↑↓ = Button for entering a value or for selecting a function.

ON/OFF = Switches the instrument on or off.

Set-up

1. Select the set-up menu by pressing **MODE**.
2. Follow the instructions on the screen to adjust language, contrast and automatic power-off timer or to erase the memory. Press **CAL** to confirm or **MODE** to exit.

The measuring electrodes should be connected to the corresponding coaxial pH/mV or EC/O₂ connectors. Automatic temperature compensation and temperature measurements are possible by plugging a Pt1000 temperature probe into the °C terminals. You can also use an electrode with built-in Pt1000. Its banana plugs should be inserted in the °C terminals. Without Pt1000, the manual temperature compensation is automatically switched on.

Important

- Never immerse a pH electrode and an oxygen electrode together in the same solution!
- A blinking decimal point warns you for unstable measurements. Wait to read the display!
- Stirring the solution during the measurements promotes the homogeneity and is obligatory!
- The instrument will refuse automatic calibration when the electrode is unstable. Insufficient stirring or a worn electrode may be the cause.

pH measurement

1. Select [pH] by pressing **MODE**. The display will immediately show the measured value according to the previous calibration. Should you want to recalibrate, press **CAL**.
2. The display shows three of the 11 buffers in memory (1.68, 2.00, 4.00, 4.01, 6.87, 7.00, 9.18, 9.21, 10.01, 12.00, 12.45). Select the proper buffers and press **CAL**. The unused buffer should be switched off.
3. Rinse the electrodes with distilled water and immerse them in the first buffer solution. Select [CALIBRATE], press **CAL** and follow the instructions on the screen until the calibration is finished.
4. After rinsing the electrodes with distilled water, immerse them in the samples and read the display.
5. Rinse the electrodes always with distilled water after use and store them in a 3...4 M KCl solution.

mV measurement

1. Select [mV] by pressing **MODE**. The display will immediately show the measured value according to the previous calibration (optional). Should you want to recalibrate, press **CAL**.
2. Immerse the electrodes in a standard solution of known potential. Select [CALIBRATE] and press **CAL**. Calibrate to the proper value and press **CAL** again or press **MODE** to reset the factory settings.
3. Rinse the electrodes always with distilled water after use and store them in a 3...4 M KCl solution.

GLP

1. Select the desired range by pressing **MODE**.
2. Press **CAL**, select GLP and press **CAL** to display a complete calibration report.

Logging data

1. Press **↓** to store the displayed value in memory. The display shows a next storage address.
2. Select the desired address and press **CAL** to enter the measured value in memory.

Recalling data

1. Press **↑** to recall the stored data. The display shows a list of stored values.
2. Select the desired address to read the stored values.
3. Press **MODE** to return to the normal measurements.

Conductivity measurement

1. Select [S/cm] by pressing **MODE**. The display will immediately show the measured value according to the previous calibration. Should you want to recalibrate, press **CAL**.
2. The display shows one of the 3 standards in memory (0.01, 0.1, 1 M KCl). Select the proper standard and press **CAL**. Also select the temperature to which all future measurements will be referred to.
3. After rinsing the cell several times with the standard solution, immerse it in that standard. The solution temperature is not so critical but should lie between 0°C and 30°C. When no Pt1000 is used, do not forget to compensate manually first! Select [CALIBRATE], press **CAL** and follow the instructions on the screen until the calibration is finished.
4. Rinse the cell several times with the sample, immerse it in that solution and read the display.
5. Rinse the cell always after use and store it in distilled water (add some detergent to keep the spongy platinum surface in perfect condition).

TDS/Salinity measurement

(C5020)

1. Select [TDS] or [SAL] by pressing **MODE**. The display will immediately show the measured value according to the previous calibration. Should you want to recalibrate, press **CAL**. Proceed as for conductivity.

Oxygen measurement

1. Select [O₂] by pressing **MODE**. The display will immediately show the measured value according to the previous calibration. Should you want to recalibrate, press **CAL**.
2. The display shows the air pressure correction and the salinity correction. Select the proper values and press **CAL**. Leave salinity correction to zero unless you are going to measure in heavily salted solutions such as e.g. sea-water (35 g/l). Select [CALIBRATE], press **CAL** and follow the instructions on the screen until the calibration is finished.
3. The electrode exposed to the air reaches an equilibrium corresponding to the partial pressure of oxygen and thus to saturation in water at the given temperature. The instrument shows the measured saturation, current, temperature, and will calibrate automatically when readings are stable.
4. After rinsing the electrodes with distilled water, immerse them in the samples and read the display. Stirring the solution during the measurements promotes the homogeneity and is **obligatory!** The advection rate must be at least 10 cm/s.
5. Rinse the electrode always after use and store it in distilled water.

%O₂ measurement

1. Select [%O₂] by pressing **MODE**. The display will immediately show the measured value according to the previous calibration. Should you want to recalibrate, press **CAL**. Proceed as for mg/l.

Temperature measurement

1. Select [°C] by pressing **MODE**. Without Pt1000, adjust the manual temperature compensation and proceed by pressing **MODE**. Should you want to recalibrate, press **CAL**.
2. Immerse the Pt1000 in a standard solution of known temperature. Select [CALIBRATE] and press **CAL**. Calibrate to the proper value and press **CAL** again or press **MODE** to reset the factory settings.

A pH electrode is active and stable only after wetting! For this purpose it must be immersed for **at least ten hours** in a 3...4 M KCl solution. During short interruptions (e.g. storage) the electrode should be immersed in a 3...4 M KCl solution. In doing this it is always kept ready for use. When the interruption is longer than a month, refill the closing cap with 3...4 M KCl and plug it on the electrode tip in order to protect the glass bulb. Before use, ensure that the reference part of the electrode is topped up with a 3...4 M KCl solution.

Avoid a low pressure inside the electrode! Therefore always remove the closure from the refilling aperture during the measurements as well as during the calibration. This allows the saltbridge solution to flow through the ceramic liquid junction and prevents contamination of the electrolyte. For the same reason, the inside level should always be higher than the outside level of the measuring solution. Close the refilling aperture again when storing the electrode.

A polluted electrode may be cleaned with a soft detergent or 0.1 M HCl. Greasy substances may be removed with acetone or alcohol (**never do this with plastic electrodes!**).

If the electrode is polluted by proteinaceous materials (such as blood), it should stand in a cleaning solution overnight and then be cleaned with distilled water before use. The pH electrode wears away by being used. If the electrode tends to respond slower and calibration becomes difficult, even after cleaning, it should be replaced by a new one.

Maintenance of metal electrodes

Metal electrodes (Pt, Ag, Au): Metal electrodes are always ready for use. During short interruptions they are immersed in distilled water. **They should be cleaned regularly:**

- Silver electrodes are immersed in a concentrated ammonia solution during one hour.
- Platinum or gold electrodes are immersed in concentrated nitric acid during one hour.

Maintenance of conductivity electrodes

A conductivity cell is active and stable only after wetting! For this purpose it must be immersed for **at least one hour** in distilled water. Rinse the cell always after use and store it in distilled water (add some detergent to keep the spongy platinum surface in perfect condition).

A polluted cell may be cleaned with a soft detergent or diluted nitric acid. Greasy substances may be removed with acetone or alcohol (**never do this with plastic electrodes!**).

Maintenance of oxygen electrodes

A dissolved oxygen electrode is active and stable only when polarised! The electrode body has almost unlimited life and requires no maintenance. After prolonged use of the electrode, it may become deactivated. An indication is that the electrode no longer responds correctly to calibration. In this case:

1. Carefully screw off the measuring head.
 2. Cautiously remove the precipitated deposits with the aid of filter paper. **Do not use grinding paper or a glass fibre!**
 3. Rinse several times with distilled water and shake carefully off the water drops.
 4. Take the measuring head, fill it with new electrolyte solution and very slowly screw it onto the electrode while holding it vertically.
 5. **Let the electrode rest for a few hours!** The electrode is now ready for use again.
- *Replace the membrane only when damaged, not when calibration is no longer possible!*

| Specifications | | C5010 | C5020 |
|------------------|--------------------------------------|--------------------------|---------------------------------|
| pH | <i>Range</i> | 0...14 pH | 0...14 pH |
| | <i>Resolution</i> | 0.01 pH | 0.01 pH |
| | <i>Accuracy</i> | 0.2% ± 1 digit | 0.2% ± 1 digit |
| | <i>Calibration</i> | 1...3 points | 1...3 points |
| | <i>Buffers</i> | 11 pre-programmed | 11 pre-programmed |
| | <i>Temperature Compensation</i> | 0...100°C | 0...100°C |
| | <i>ISO-pH</i> | 6...8 pH | 6...8 pH |
| | <i>Slope</i> | 80...120% | 80...120% |
| mV | <i>Range</i> | ±1000 mV | ±1000 mV |
| | <i>Resolution</i> | 1 mV | 1 mV |
| | <i>Accuracy</i> | 0.2% ± 1 digit | 0.2% ± 1 digit |
| | <i>Calibration</i> | 1 point | 1 point |
| CONDUCTIVITY | <i>Range</i> | 0...100 mS/cm | 0...1000 mS/cm |
| | <i>Resolution</i> | 0.01 µS/cm | 0.01 µS/cm |
| | <i>Accuracy</i> | 1% f.s. of range | 1% f.s. of range |
| | <i>Calibration</i> | 1 point | 1 point |
| | <i>Standards</i> | 0.01/0.1M KCl | 0.01/0.1/1M KCl |
| | <i>Cell Constant</i> | 1 cm ⁻¹ ± 30% | 0.1/1/10 cm ⁻¹ ± 30% |
| | <i>Temperature Compensation</i> | 0...100°C | 0...100°C |
| | <i>Reference Temperature</i> | 20° or 25°C | 20° or 25°C |
| | <i>Temperature Coefficient</i> | natural waters (EN27888) | natural waters (EN27888) |
| SALINITY | <i>Range</i> | | 0...70 ppt |
| | <i>Resolution</i> | | 0.1 ppt |
| TDS | <i>Range</i> | | 0...100 g/l |
| | <i>Resolution</i> | | 0.1 mg/l |
| DISSOLVED OXYGEN | <i>Range</i> | 0...20 mg/l (0...200%) | 0...20 mg/l (0...200%) |
| | <i>Resolution</i> | 0.05 mg/l (0.5%) | 0.05 mg/l (0.5%) |
| | <i>Accuracy</i> | 1% ± 1 digit | 1% ± 1 digit |
| | <i>Calibration</i> | 1 point | 1 point |
| | <i>Temperature Compensation</i> | 0...50°C | 0...50°C |
| | <i>Salinity Compensation</i> | 0...40 ppt | 0...40 ppt |
| | <i>Air Pressure Compensation</i> | 800...1200 hPa | 800...1200 hPa |
| TEMPERATURE | <i>Range</i> | 0...100°C | 0...100°C |
| | <i>Resolution</i> | 0.1°C | 0.1°C |
| | <i>Accuracy</i> | 0.5°C | 0.5°C |
| | <i>Calibration</i> | 1 point | 1 point |
| INPUTS | <i>pH/mV</i> | BNC, 10 ¹² Ω | BNC, 10 ¹² Ω |
| | <i>Conductivity/Dissolved Oxygen</i> | BNC | BNC |
| | <i>Temperature</i> | 2 banana, for Pt1000 | 2 banana, for Pt1000 |
| STORAGE MEMORY | <i>Data sets</i> | | 300 |
| DISPLAY | <i>LCD</i> | 122x32 pixels | 122x32 pixels |
| | <i>White Backlight</i> | ✓ | ✓ |
| | <i>Languages</i> | EN, NL, FR, DE | EN, NL, FR, DE |
| POWER SUPPLY | <i>Mains</i> | 210...250 VAC, 50/60 Hz | 210...250 VAC, 50/60 Hz |
| | <i>Low Voltage</i> | 9...15 VDC | 9...15 VDC |
| | <i>Batteries</i> | 4x1.2 V, NiMH | 4x1.2 V, NiMH |
| DIMENSIONS | <i>WxDxH</i> | 10x20x4 cm | 10x20x4 cm |
| WEIGHT | <i>Meter</i> | 350 g | 350 g |

| Temp. | 1.68 | 4.01 | 6.87 | 9.18 | 12.45 | <i>NIST buffers</i> | |
|-------|-------|-------|-------|-------|--------|---------------------|--|
| 0° C | 1.666 | 4.003 | 6.984 | 9.464 | 13.423 | (pH) | |
| 5° C | 1.668 | 3.999 | 6.951 | 9.395 | 13.207 | | |
| 10° C | 1.670 | 3.998 | 6.923 | 9.332 | 13.003 | | |
| 15° C | 1.672 | 3.999 | 6.900 | 9.276 | 12.810 | | |
| 20° C | 1.675 | 4.002 | 6.881 | 9.225 | 12.627 | | |
| 25° C | 1.679 | 4.008 | 6.865 | 9.180 | 12.454 | | |
| 30° C | 1.683 | 4.015 | 6.853 | 9.139 | 12.289 | | |
| 40° C | 1.694 | 4.035 | 6.838 | 9.068 | 11.984 | | |
| 50° C | 1.707 | 4.060 | 6.833 | 9.011 | 11.705 | | |
| 60° C | 1.723 | 4.091 | 6.836 | 8.962 | 11.449 | | |
| 70° C | 1.743 | 4.126 | 6.845 | 8.921 | | | |
| 80° C | 1.766 | 4.164 | 6.859 | 8.885 | | | |

| Temp. | 2.00 | 4.00 | 7.00 | 9.21 | 10.01 | 12.00 | <i>Standard buffers</i> | |
|-------|------|------|------|------|-------|-------|-------------------------|--|
| 5° C | 1.85 | 3.99 | 7.08 | 9.45 | 10.24 | 12.56 | (pH) | |
| 10° C | 1.89 | 3.99 | 7.06 | 9.38 | 10.18 | 12.35 | | |
| 20° C | 1.96 | 3.99 | 7.02 | 9.26 | 10.06 | 12.14 | | |
| 25° C | 2.00 | 4.00 | 7.00 | 9.21 | 10.01 | 12.00 | | |
| 30° C | 2.03 | 4.01 | 6.99 | 9.16 | 9.87 | 11.89 | | |
| 40° C | 2.08 | 4.03 | 6.98 | 9.06 | 9.89 | 11.71 | | |
| 50° C | 2.12 | 4.06 | 6.97 | 8.99 | 9.83 | 11.52 | | |
| 60° C | 2.14 | 4.09 | 6.98 | 8.93 | 9.79 | 11.33 | | |

| Temp. | 0.01 M | 0.1 M | 1 M | <i>Conductivity standards</i> | | | |
|-------|--------|-------|-------|-------------------------------|--|--|--|
| 0° C | 0.776 | 7.15 | 65.4 | (mS/cm) | | | |
| 5° C | 0.896 | 8.22 | 74.1 | | | | |
| 10° C | 1.020 | 9.33 | 83.2 | | | | |
| 15° C | 1.147 | 10.48 | 92.5 | | | | |
| 20° C | 1.278 | 11.67 | 102.1 | | | | |
| 25° C | 1.413 | 12.88 | 111.8 | | | | |
| 30° C | 1.552 | 14.12 | 121.7 | | | | |

| Temp. | 124 | 358 | <i>Redox standards</i> | | | |
|-------|-----|-----|------------------------|--|--|--|
| 0° C | 156 | 388 | (mV) | | | |
| 5° C | 148 | 381 | | | | |
| 10° C | 141 | 376 | | | | |
| 15° C | 135 | 370 | | | | |
| 20° C | 130 | 364 | | | | |
| 25° C | 124 | 358 | | | | |
| 30° C | 122 | 352 | | | | |
| 35° C | 116 | 346 | | | | |
| 40° C | 99 | 340 | | | | |



DECLARATION OF CONFORMITY

We declare under our sole responsibility that the product

Multi-channel analysers
content of the type numbers
C5010, C5020

to which this declaration relates is in conformity

with the following standards

EN61010

LOW VOLTAGE DIRECTIVE 73/23/EEG

EN50081-1

EN50082-1

EN60555-2

EMC DIRECTIVE 89/336/EEG

Turnhout, February 02, 2011

A handwritten signature in black ink, appearing to be 'L. G.', is located below the date.

on behalf of Consort

WARRANTY CERTIFICATE

This instrument (excluding all accessories) is warranted against defective material and workmanship for a period of thirty-six (36) months from the date of shipment ex factory. Accessories and breakable items such as electrodes are not warranted unless proven to be defective before shipment.

The original purchase order numbers, Consort invoice numbers and serial numbers of the products must be provided.

CONSORT will repair all defective equipment returned to it during the warranty period without charge (CIF Turnhout prepaid by sender), provided the equipment has been used under normal laboratory conditions and in accordance with the operating limitations and maintenance procedures described in the instruction manual and when not having been subject to accident, alteration, misuse or abuse.

If the products have been used with or have come into contact with fluids, an MSDS (material safety data sheet) must be supplied prior to issuing a return authorisation.

A return authorisation must be obtained from **CONSORT** before returning any product for warranty repair on a freight prepaid basis!

Consort

www.consort.be

Parklaan 36 • B-2300 Turnhout • Belgium

Tel: (+32) (0)14 41 12 79

Fax: (+32) (0)14 42 91 79

Sales: sales@consort.be

Support: support@consort.be

Information: info@consort.be